


**KeriFlex®**




**KeriFlex** Finger Joint Implants from KeriMedical are flexible, one-piece, hinged implants manufactured from silicone designed to replace damaged metacarpophalangeal and proximal interphalangeal joints (MCP and PIP, respectively). They provide you with a strong and durable option when treating painful MCP and PIP joints in your patients.

## 12 Sizes Available

Hinge dimensions shown below

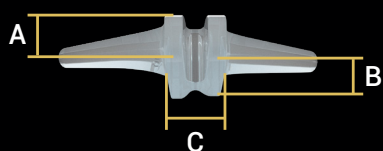

PIP - 5 sizes of implants			
Size	A	B	C
1	9 mm	6.5 mm	6 mm
2	10 mm	7 mm	6.9 mm
3	10.8 mm	7.8 mm	7.6 mm
4	12 mm	8.5 mm	8.2 mm
5	13 mm	9.4 mm	9.1 mm


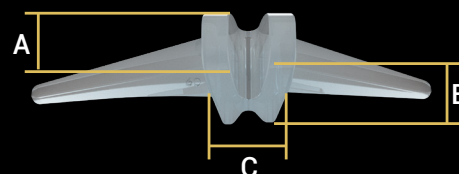


MCP - 7 sizes of implants			
Size	A	B	C
00	9 mm	8.7 mm	6.7 mm
10	10 mm	9.4 mm	7.5 mm
20	10.8 mm	10.2 mm	8 mm
30	12 mm	11.4 mm	9 mm
40	13 mm	12.4 mm	10 mm
50	13.8 mm	13.5 mm	10.5 mm
60	15 mm	14.6 mm	11.6 mm

## Anatomically Neutral Implant Design

- Pre-flexed, joint-hinge design mimics how our hands naturally rest (15° PIP/ 30° MCP), while decreasing stress on the implant
- Provides flexion in the hinge region
- Square shape hinge is reduced in height to limit conflicts with extensors apparatus and increased in width to improve lateral stability
- Anatomic shape offers a better fit
- 12 total sizes (5 PIP/ 7 MCP) to accommodate varying patient anatomies
- Midsection maintains proper joint space and alignment with good lateral stability with no flexion/extension restriction

Full flexion/extension (0° to 90°)

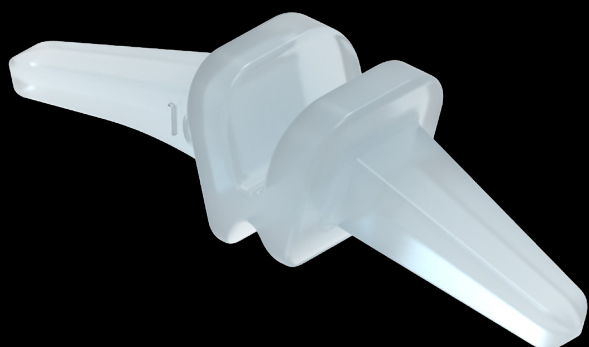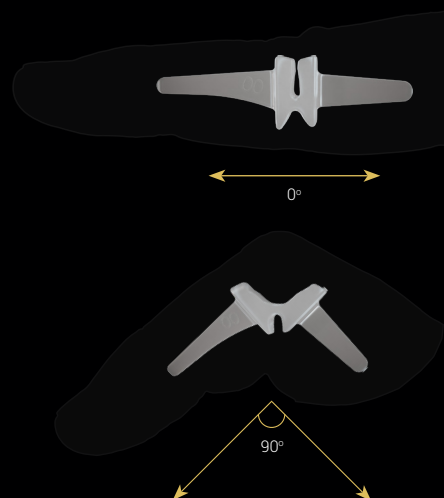


## Durable Material. Enhanced Experience.

- Superior quality silicone which performed verified resistance testing up to 10 million flexion/extension cycles with no rupture<sup>1</sup>
- Short procedure
- Easily performed in an ASC setting
- Can be performed under local anesthesia

## Simplified Instrumentation

- Easy-to-use, ergonomically designed instruments
- A compact PIP & MCP set offers everything in one
- 2 progressive rasps per indication (4 total)

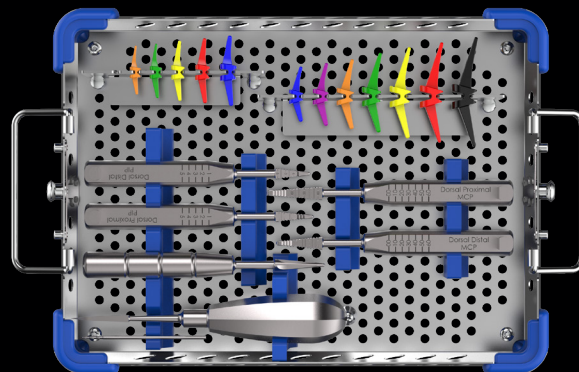

Validated sterilization cycle 132°C (270°F) / 4 minutes is also validated according to ANSI/AAMI ST79 and EN ISO 17665

<sup>1</sup>. Data on file with KeriMedical

**INDICATIONS:** The KeriFlex® MCP and KeriFlex® PIP Finger Joint Prostheses are indicated for cementless replacement of the metacarpophalangeal (MCP) and interphalangeal (PIP) joints, respectively, where disabled by rheumatoid, degenerative or traumatic arthritis.

### Legal Manufacturer:

KeriMedical SA  
Route des Acacias, 45a  
1227 Geneva - Switzerland  
Tel: +41 58 255 01 30  
Fax: +41 58 255 01 40  
customerservice@kerimedical.com

### Distributed by:

This material is intended for health care professionals. Distribution to any other recipient is prohibited. For product information, including indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert.

KeriFlex is a trademark of KeriMedical SA

©2024 Medartis Inc.

KERI-0014 V0 04/2024