

medartis®

PRODUCT INFORMATION

All-in-One Staple



APTUS®
Foot



Ready. Steady. Staple.
The Swiss Pocket Staple.


Ready. Steady. Staple.

The pocket size sterile staple kit with implant and all implant-specific instruments in a lightweight packaging of 19 g and 66 x 99 x 11 mm.

A further step towards comprehensive hallux treatment and bone fixation in the foot.


The Swiss Pocket Staple

- Complete
- Functional
- Sterile and compact on the shelf

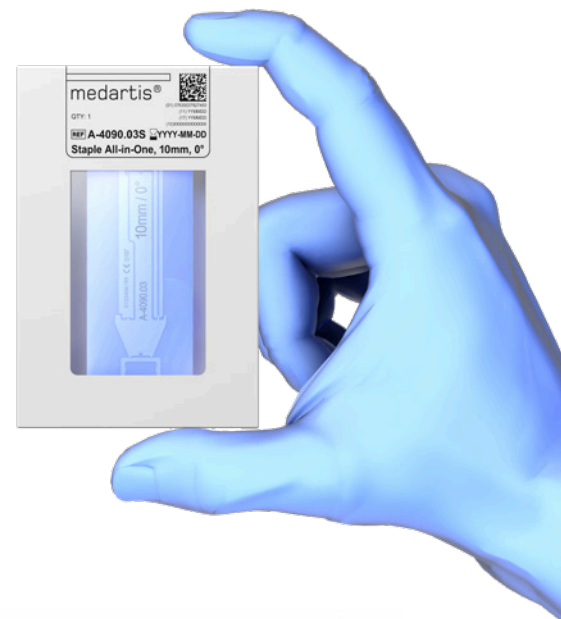

Implant and Instruments

- 1-piece precision laser cut from stainless steel sheet material
- Implant and instruments can be easily separated
- Different staple designs to meet anatomical needs and surgeon's preference
- Swiss engineering and manufacturing

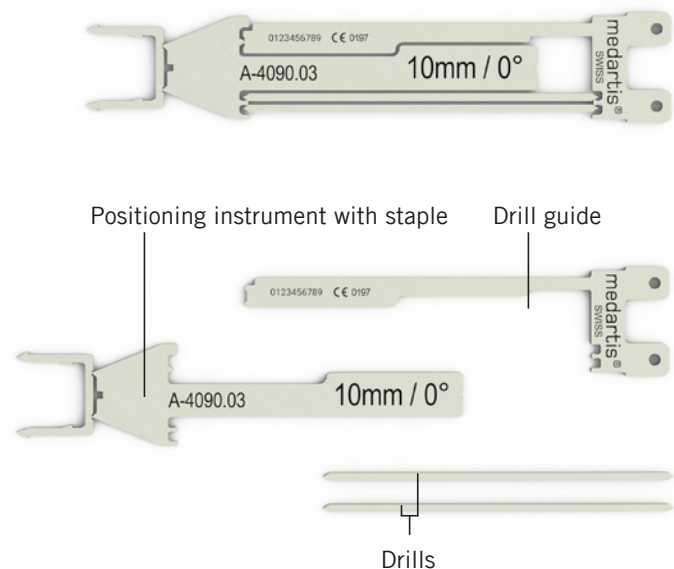

Scale 1:1

Specifications

- Available in two widths (8 mm and 10 mm) and two angles (0° and 26°)
- Leg length 10 mm
- Barbed legs can prevent migration of the staple
- Staple design and drill guide allow for compression
- Two drills and the corresponding drill guide are included
- Positioning instrument for precise insertion

Values

- Saving reprocessing capacity and resources
- Lean material and supply management – 1 item, 1 price
- No additional costs for consumables such as K-wires for drilling
- Reducing OR idle time



Staple All-in-One
8 mm, 0°



Staple All-in-One
8 mm, 26°



Staple All-in-One
10 mm, 0°



Staple All-in-One
10 mm, 26°

Ordering Information

Material: Stainless steel (ASTM F139)
Staple thickness: 1.0 mm

Staples All-in-One



A-4090.01S



A-4090.03S



A-4090.02S



A-4090.04S

Art. No.		Description (Width, Angle)	Pieces/Pkg
A-4090.01S	STERILE	08 mm, 0°	1
A-4090.02S	STERILE	08 mm, 26°	1
A-4090.03S	STERILE	10 mm, 0°	1
A-4090.04S	STERILE	10 mm, 26°	1

For further information regarding APTUS Foot visit www.medartis.com.

Scale 1:1

FOOT_07000006_v0 / © 2019-09, Medartis AG, Switzerland. All technical data subject to alteration.

MANUFACTURER & HEADQUARTERS

Medartis AG | Hochbergerstrasse 60E | 4057 Basel/Switzerland
P +41 61 633 34 34 | F +41 61 633 34 00 | www.medartis.com

USA

Medartis Inc. | 224 Valley Creek Boulevard, Suite 100 | Exton, PA 19341
P +1 610 961 6101 | Toll free 877 406 BONE (2663) | F +1 610 644 2200

SUBSIDIARIES

Australia | Austria | Brazil | France | Germany | Japan | Mexico | New Zealand | Poland | UK | USA

For detailed information regarding our subsidiaries and distributors, please visit www.medartis.com



Disclaimer: This information is intended to demonstrate the Medartis portfolio of medical devices. A surgeon must always rely on her or his own professional clinical judgement when deciding whether to use a particular product when treating a particular patient. Medartis is not giving any medical advice. The devices may not be available in all countries due to registration and/or medical practices. For further questions, please contact your Medartis representative (www.medartis.com). This information contains CE-marked products.
For US only: Federal law restricts this device to sale by or on the order of a physician.